

# Urgent Stabilisation for Children's Femoral Fractures



# Transfer of an injured child



- Clinically evident fractured SOF
- Analgesia
- Traction Splint
  - Kendrick Traction Device maintained by skin traction
  - Thomas splint
- Xray ?
- Trauma Blue Transfer via NWAS


## Benifits

- Manage pain
- Limit blood loss
- Mitigate Neurovascular injury

# Equipment

- Adhesive Skin traction
- KTD
- Gamgee padding/ wound pads
- Bandage
- Tongue depressors

# Adhesive Skin Traction

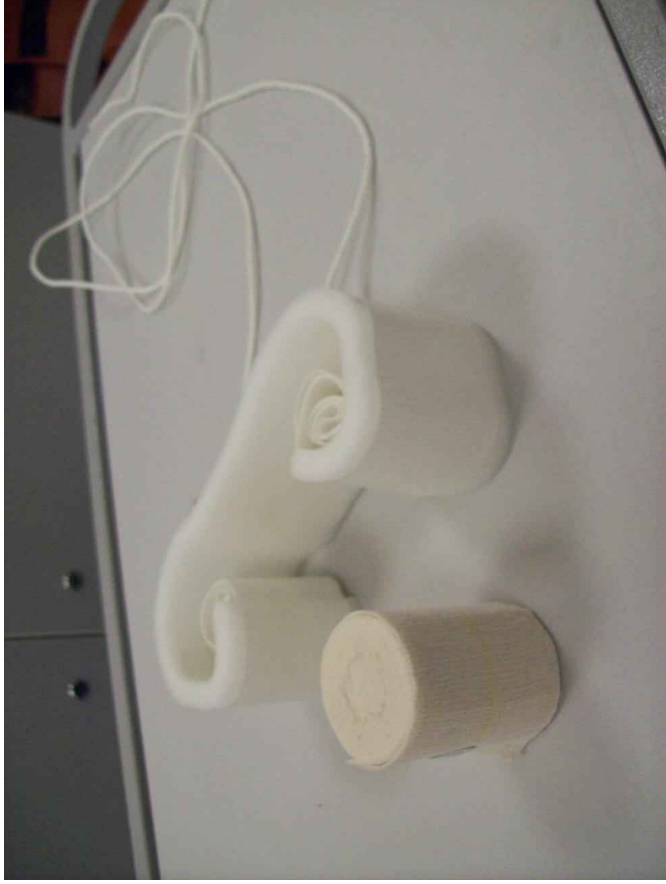

# KTD in children?



# Components of KTD

- Use the KTD pole and groin strap
- Supportive straps can be used
- Don't use the ankle strap

Swollen,  
Shortened and Internally rotated






# Neurovascular Status



# Analgesia

- Opiate IN/ IV
- Entonox (demand or free-flow)
- Ketamine
- Femoral nerve block
- RSI (for other injuries)

## KTD groin strap

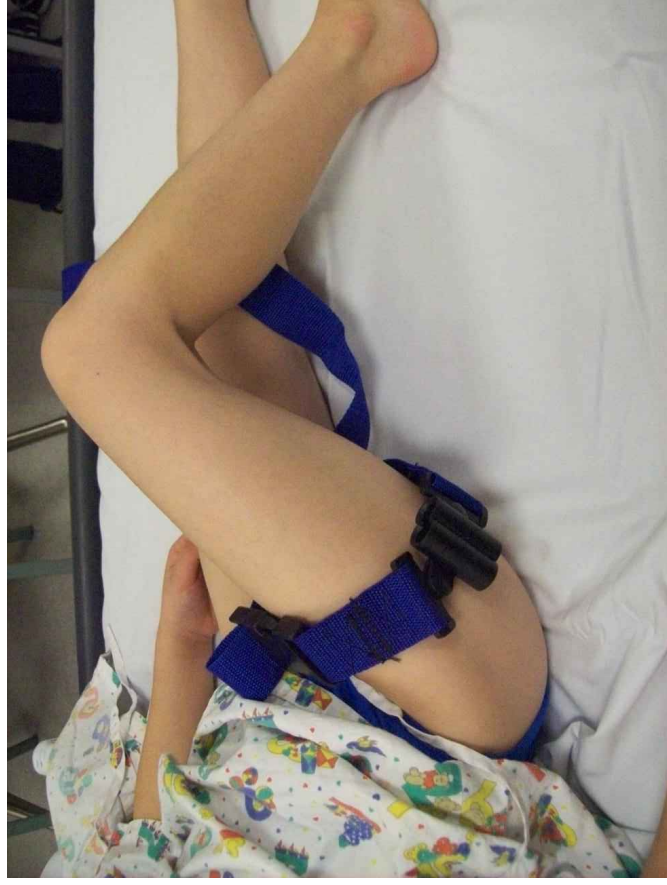

After analgesia apply traction



# Heel pad



# Apply adhesive straps



# Bandage






## Apply traction via cords





## Insert poles in the groin strap

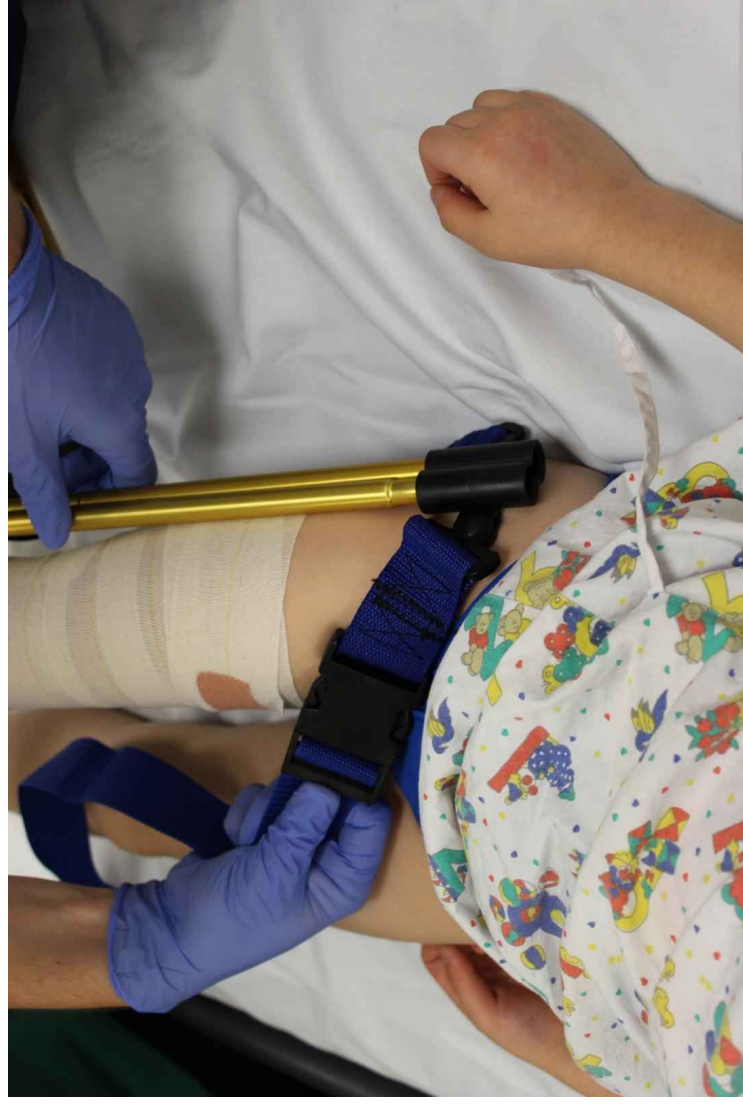


Push and secure groin strap



# Attach to end bar



## Secure with Figure of '8' knot

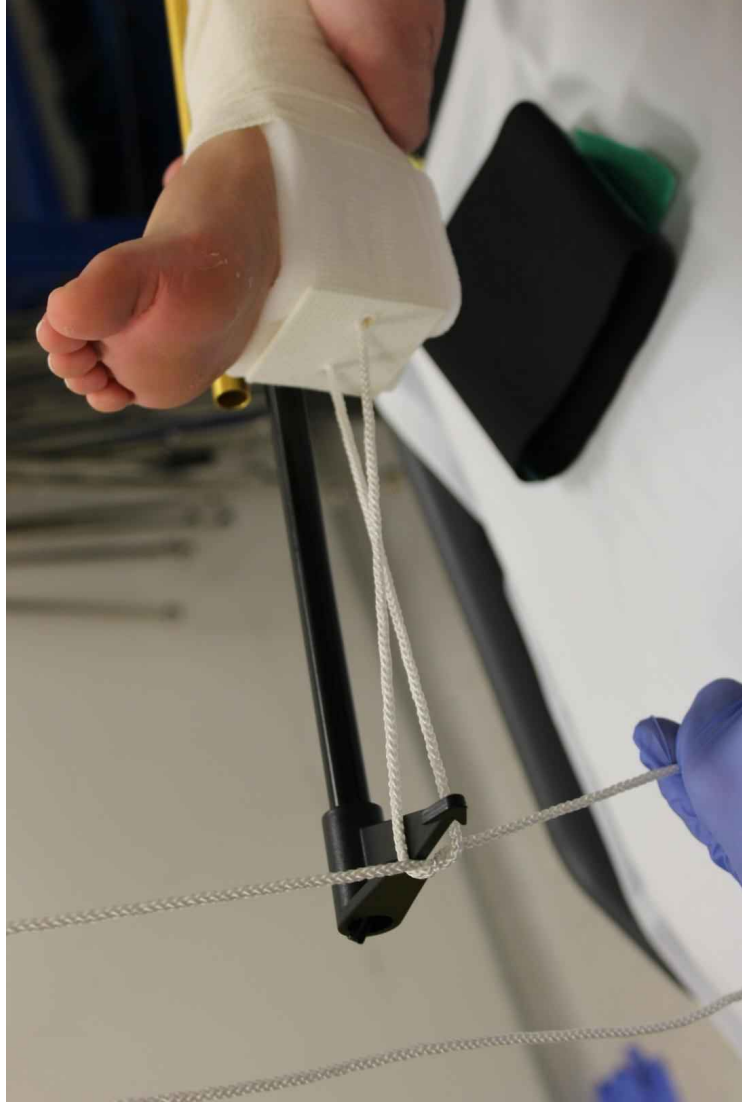

# Pad with Gamgee for Transfer



# Bandage splint





# Spanish Windlass



?





## Conclusions

- Simple, Quick and Effective
- X Ray not needed in obvious deformity
- Skin traction transfers to Thomas Splint
- Availability NWAS/ TU